

Sample Cutter

Tall Universal Tray



1" Template – **P/N 90-07-T**

15 mm Template – **P/N 90-07-TM**

Max Tray Depth: 5"

Specifications:

- Stainless Steel template allows for error-free tray sample cutting.
- Adjustable height for trays up to 5" (130mm) deep.
- Dimensions: 13"L x 6.5"W x 10.5"H (330mm x 165mm x 270mm)
- Instrument weight: 10 lbs. (4.5 kg)

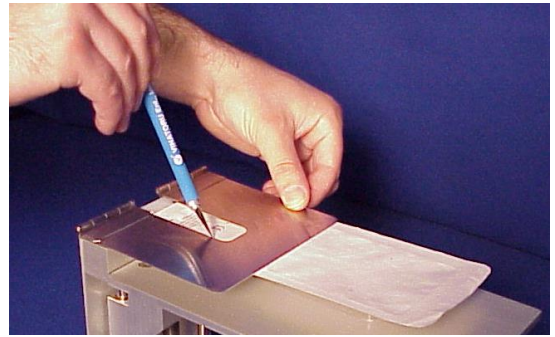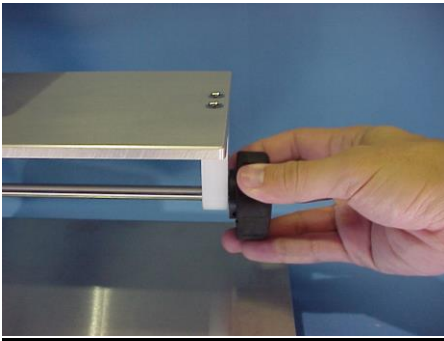

1. Adjust the tray support height on both the universal sample cutter and the HST 180T instrument to suit the tray to be tested.

2. Determine location on the tray from where the sample will be cut. Place the tray flange under the hinged template and push down trapping the seal over plastic insert. Using the supplied trimming knife, cut out the sample sides. Make sure the blade used is sharp so that the lid is cut during the first attempt.

For best results use:

- Trimming Knife **P/N 90-00-01** - \$5
- Replacement Blades (100 pcs.) **P/N 90-00-00** - \$25
- USC II - Plastic Insert **P/N 10-06-141** - \$7.50

Price

Universal Tray Sample Cutter Tall **P/N 90-07-T** or **P/N 90-07-TM** - \$ 740

Note: Certificate of Compliance available - \$20

If needed, please specify at time of purchase.

May be purchased online at www.vinatoru.com.